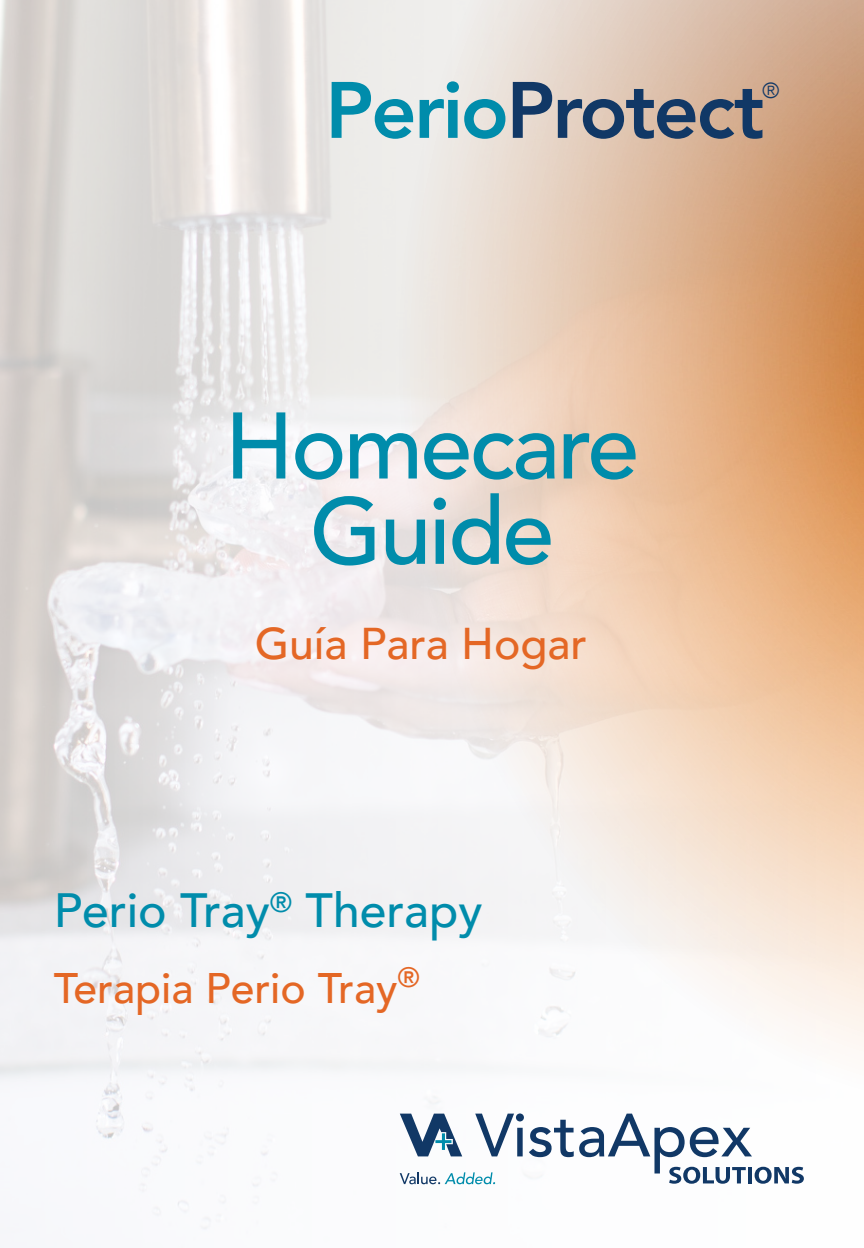
PerioProtect®

Homecare Guide

Guía Para Hogar

Perio Tray® Therapy

Terapia Perio Tray®

 **VistaApex**
SOLUTIONS
Value. Added.

Instructions for Use

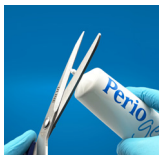
Figure A

Figure B

Figure C

Visit the
PerioProtect®
Resource Center
(See QR code
on back)

1. Cut the tip of your Perio Gel® tube [A].
2. Hold the Perio Tray® open side up [B].
3. Apply a small, thin line across the tooth indentations [C].
4. Place the tray into your mouth so the tray covers all the teeth and fully seats in place. Spit out any excess medication.
5. Wear your trays for the time prescribed by your doctor. (Typically, 15 minutes per session)
6. After use, remove the trays. Brush your teeth to remove any excess medication.
7. Use your toothbrush to clean the trays and rinse them thoroughly with cold water. **Do not use hot water.**
8. Gently towel-dry your trays and return them to the storage case.

Perio Tray® Wearing Instructions

(Clinician to complete)

For patients diagnosed with periodontitis

- ☐ 3 times a day for 15 minutes
- ☐ 2 times a day for 15 minutes
- ☐ 1 time a day for 15 minutes
- ☐ _____

For patients diagnosed with gingivitis and for maintenance

- ☐ 2 times a day for 15 minutes
- ☐ 1 time a day for 15 minutes
- ☐ _____

It is very important that you continue to brush and floss your teeth. If you experience any adverse reactions during Perio Tray® Therapy, discontinue use and contact your dentist.

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PerioProtect, LLC